

Magnetic Switch

BE9023-P02



Features

Flexible
Works for doors and windows

Easy to set up
Pre-mounted adhesive tape

Visit compatible
Connects to the Visit system


Monitors your door and window.

The Magnetic Switch from Bellman & Symfon is specially designed for the Visit System and alerts you when a door or window is opened.

Just mount the magnetic switch on the door or window frame and connect it to the Visit telephone transmitter. When the magnets are separated, the transmitter signals your Visit receiver that alerts you with sound, flashes or vibrations (depending on the receiver).

Technical specifications

- Dimensions 25 x 62 x 13 mm 1" x 2.5" x 0.5"
- Weight 25 g, 0.9 oz.
- Connector 3.5 mm mono plug
- Cable length 0.5 m, 1.6'
- Installation Make sure the distance between the magnets is no more than 1cm when the door is closed.
- Colour White
- For indoor use only

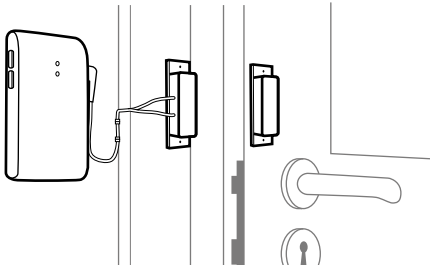

Approvals

MD This product is a medical device in conformity with Regulation (EU) 2017/745 (Medical Device Regulation).

© 2025, Bellman & Symfon AB. All rights reserved.
Bellman, Bellman & Symfon and the Bellman & Symfon logo are registered trademarks of Bellman & Symfon AB. This leaflet is subject to change without prior notice.

BE9023-P02_003DSH1.0_EN bellman.com