

Mobile Phone Transceiver

BE1433



Connects Visit to your mobile phone, tablet or smartwatch.

The Mobile Phone Transceiver makes it easier than ever to keep tabs on your home. It connects the Visit system to your mobile phone, tablet, or smartwatch via Bluetooth. Once you have installed the Visit app, you will start to receive Visit notifications on your mobile device. In addition, your Visit receiver will alert you to incoming calls and messages on your mobile device.

It complements your existing receiver by providing additional information on your alerts, allowing you to customize your notifications, browse the notifications history and get an overview of your connected devices.

The Mobile Phone Transceiver can also alert to your landline telephone and has a trigger port that can be used for a contact mat, doorbell, or magnetic switch accessory.



Features

Bluetooth powered
Up to 40 m coverage, open field


App controlled
Available for iOS & Android


Sends Visit alerts to mobile devices
Displays Visit alerts on your phone, tablet or smartwatch

Alerts for mobile calls and messages
Signals your Visit receiver when you get incoming calls or messages

MD

CE-marked medical device under EU MDR 2017/745

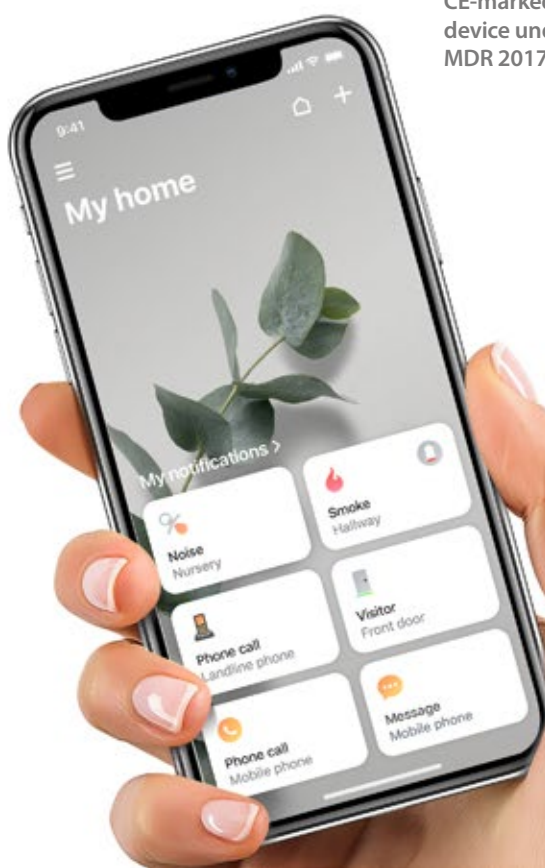

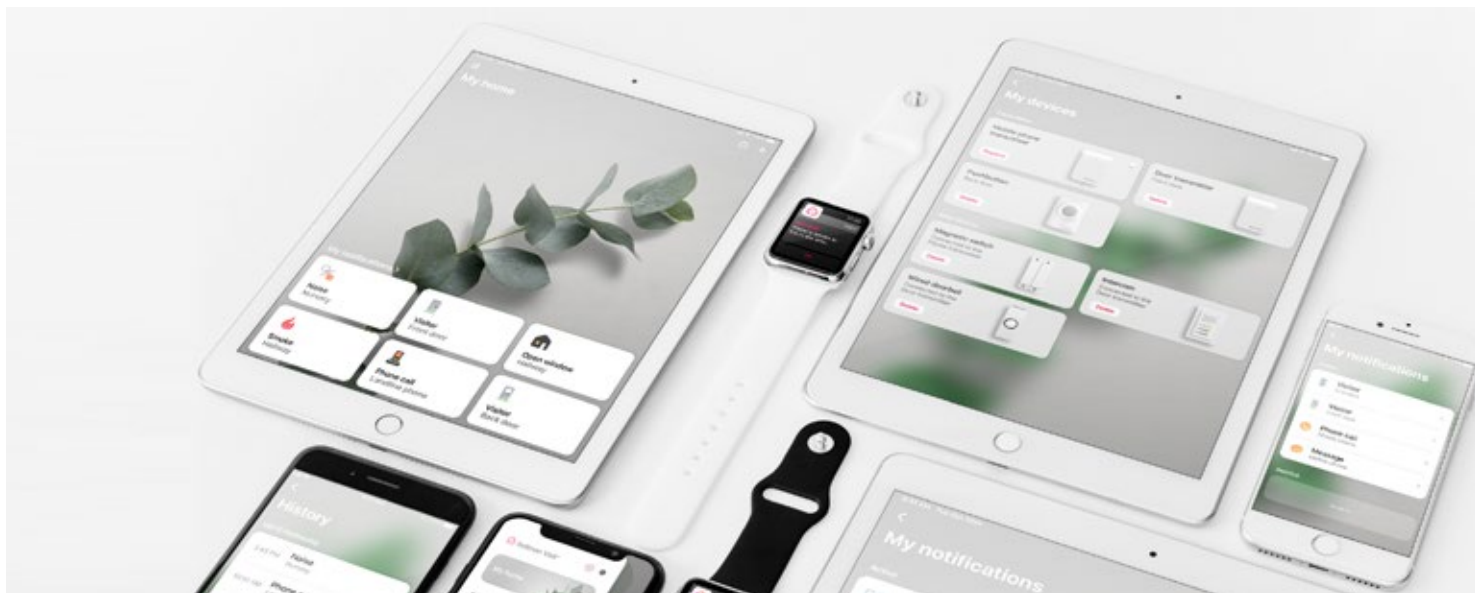


Size and weight

Height:	100 mm, 4.0"
Width:	65 mm, 2.6"
Depth:	27 mm, 1.1"
Weight:	198 g, 6.98 oz. incl. batteries



Buttons and connectors

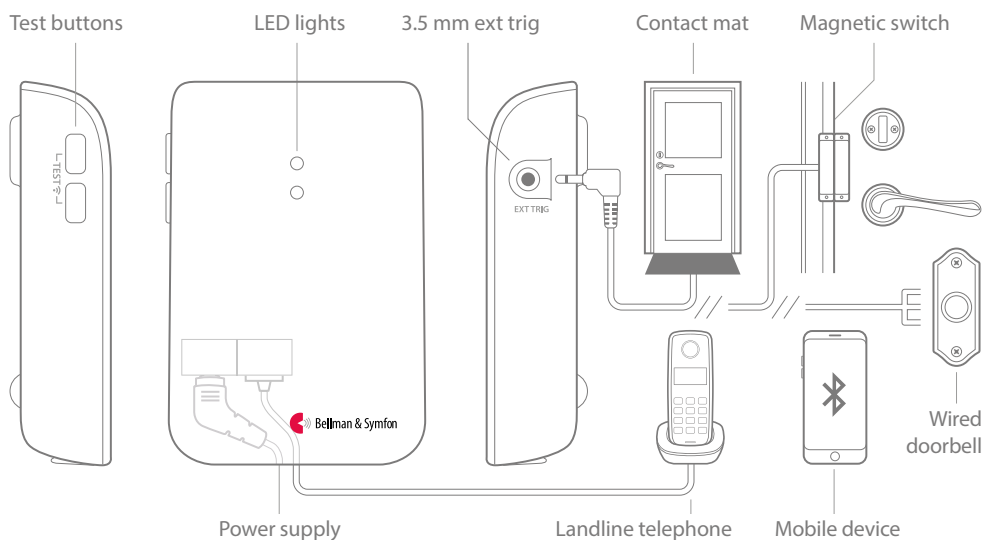

Connectivity

- RJ11 input to analogue telephone network
- 3.5 mm ext. trig port for accessories
- Power supply jack
- Bluetooth 5.0

Activation

- A mobile device via Bluetooth
- The test buttons and the landline telephone
- The contact mat or magnetic switch accessory
- A wired doorbell connected to the ext. trig. port

Compatibility

- iOS15 and later | Android 8 and later
- iPhone 6S and later | Modern Android devices

Environmental requirements

- For indoor use only.
- Temperature: 0° to 35° C, 32° to 95° F
- Relative humidity: 15% – 90% (non condensing)

Power and battery

- Mains power: 5 V DC / 1000 mA
- Backup batteries: 2 x 1.5 V AA Alkaline

Frequency and coverage

- ISM frequency: 868.30 MHz
- ISM coverage: Up to 250 m (275 yd)
- Bluetooth frequency: 2400 – 2483.5 MHz
- Bluetooth coverage: Up to 40 m

Number of radio keys

- 64 as standard.

Maintenance and cleaning

- Use a soft, lint-free cloth. Avoid getting moisture in openings. Don't use household cleaners.
- Not reusable.

Approvals



This product is a medical device in conformity with Regulation (EU) 2017/745 (Medical Device Regulation).

In the box

- BE1433 Mobile Phone Transceiver
- Power adapter and power cable
- Pre-mounted telephone cable
- 2 x 1.5V AA alkaline batteries
- Telephone splitter
- Screw and plug



Accessories

The following accessories are available:

- BE9023-P02 Magnetic switch
- BE9026-P02 Contact mat
- BE9253 External trigger cable (3.5 mm connector)

© 2025 Bellman & Symfon AB. All rights reserved

Bellman, Bellman & Symfon and the Bellman & Symfon logo are registered trademarks of Bellman & Symfon AB. This leaflet is subject to change without prior notice.

BE1433_053DSH1.0_EN

bellman.com