

# Amicus pager

Model: BE2350



## Features

The Amicus pager is robust, lightweight and easy to use. It operates for 24 hours on a single charge and has a customizable inlay. In addition, the desktop charger features battery backup and an outlet for a bed shaker.



### Clear notifications

Alerts the caregiver with a combination of powerful vibrations, colored LEDs and an optional sound signal.



### Auto configuration

Handles up to four Amicus transceivers and automatically adapts the signal patterns to the system configuration.



### Nighttime alerts

Comes with a desktop charger that is equipped with an outlet for a bed shaker to enable nighttime alerts.



### Smart tracking

Immediately warns for system faults or if the pager moves out of range, providing piece of mind for the caregiver.



### Long lasting power

Operates for a straight 24 h even during the busiest work shift and alerts for incoming alarms during charging.



### Easy to use

Features clean and functional Scandinavian design with great ergonomics and a user friendly interface.

# Technical specifications



## Pager dimensions and weight

Height:	91 mm
Width:	55 mm
Depth:	23 mm
Weight:	60 g incl. battery





## Charger dimensions and weight

Height:	110 mm
Width:	93 mm
Depth:	110 mm
Weight:	195 g incl. battery backup



## Pager and charger overview

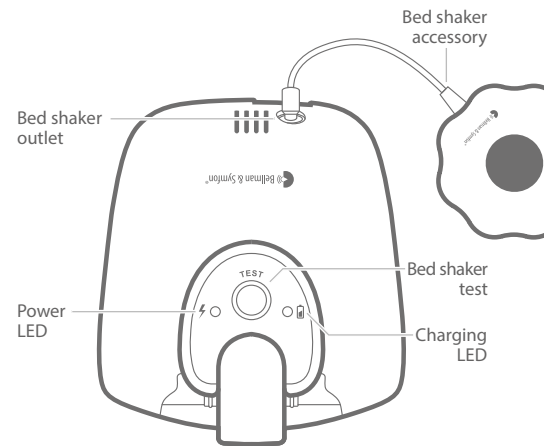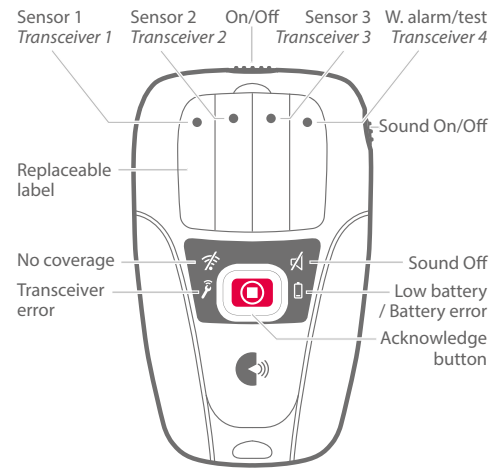

## Pager power and battery

- Battery power**  
 1 x 1.2 V AAA rechargeable NiMH battery. Use the supplied rechargeable battery model HFR-AAA900 900 mAh by Shenzhen Highpower Technology Co. Ltd.
- Operating time**  
 24 hours
- Power consumption**  
 20 mA
- Charging time**  
 Up to 6 hours

## Charger power and battery

- Mains power**  
 Input: 100-240V~50/60Hz 0.5A,  
 Output: 7.5V = 1600mA.  
 Use the supplied adaptor model SFP0751600P made by Helms-Man.
- Battery power**  
 4 x 1.2 V AAA rechargeable NiMH batteries.  
 Use the supplied rechargeable batteries model 80AAAHC 800 mAh by GP.
- Power consumption**  
 200 mA

## Frequency and coverage

- Radio frequency: 869.2375 MHz
- Coverage: Up to 600 m line of sight. The range is reduced by walls and large objects and can be affected by other radio transmitters.

## Charger connections

- 7.5 V power cord
- 3.5 mm tele jack for bed shaker  
 Features a bed shaker test button located on the top

## Environment and regulations

- For indoor use only
- Temperature: 59° to 95° F (15°-35° C)
- Relative humidity: 5% -95%, non condensing
- Complies with CE, WEEE, R&TTE and RoHS

## In the box

- BE2350 Amicus pager
- BE1251 Pager charger
- Mains adaptor
- User manual



## Accessories

The following accessory is available:

- BE1270 Bed shaker for nighttime alerts

© and TM 2015, Bellman & Symfon AB

Amicus, Bellman and Bellman & Symfon are registered trademarks owned by Bellman & Symfon AB, all rights reserved. All product designs are patented. This leaflet is subject to change without prior notice.

BE2350\_003DSH001

[bellman.com](http://bellman.com)